



European Quality Parts

Work Mode Chart

	134° C Package	134° C Non Package	121° C Package	121° C Non Package	134B+
Sterilizing Temperature	134	134	121	121	134
Sterilizing Pressure	2.10	2.10	1.10	1.10	2.10
Sterilizing Time	6	4	20	16	18
Drying time min	15	15	15	15	15
Quick Drying time	1/5	1/5	1/5	1/5	1/5
Total time min	43 (29/33)	35 (21/25)	57 (43/47)	47 (33/37)	55 (41/45)

FENG 23L / SEA 23L-B-LED

- + B type fully automatic front-loading autoclave powered by triple vacuum cycle
- + Temperature Selection - 121° C / 134° C/ B+
- + Double door lock system for safety
- + 5 Programs for sterilization of all types of instruments - Solid/Hollow, Porous/Non-Porous etc.
- + Reduced full cycle time
- + Bow and Dick test qualified
- + User-friendly Control Panel
- + 2 Storage tanks for distilled and waste water
- + Data storage facility with USB functioning (2000 cycles)
- + Steam generator cleaning program to increase the efficiency of the autoclave
- + Chamber size 250 mm x 450 mm
- + High vacuum pump level -92 kPa
- + Optional micro printer

23 Ltr B-Type front loading Fully automatic with USB RECORDER




Aluminum Foil Installation



High Vacuum Pump



Stainless Steel Connectors

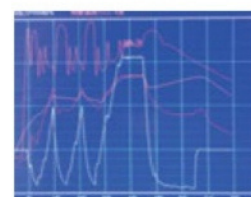

Reduced Full Cycle Time

Rahul Sali

Phone no. +91 96990 89894
E-mail: rahul.sali@lifechemsystems.com

Sourav Sali

Phone no. +91 88799 23116
E-mail: sourav.sali@lifechemsystems.com