



*Imaging  
Healthcare  
Signage*

# VMS MediPro Solution



# VMS MediPro Solution


- The **VMS Digital Radiology Solution** offers a holistic, end-to-end management option
- Ideally-suited for Hospitals, Diagnostic centres with MRI, CT-Scan, Speciality clinics for Orthopaedic, Cardiology, Dental, Endoscopic, Sonography, etc. looking to adopt new digital imaging technologies or radiology centres with few per-day exposures
- Provides easy integration to all digital modalities (CT scan, MRI, Ultrasound, Digital X-Rays, Nuclear Medicine)
- Makes the most of the digital workflow while maintaining a perfect balance between image quality and convenience
- Provides significant costs reduction as well as helps minimise the environmental impact on a daily basis.



Manage



Print

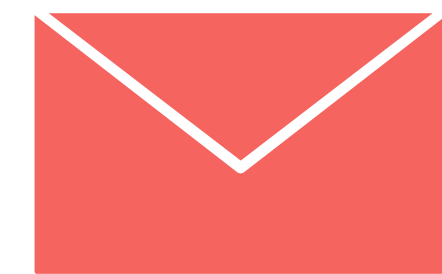

Share



Manage

VMS

MediPro™

The Most Intuitive DICOM Software



# VMS MediPro DICOM Software

MANAGE


- **VMS MediPro** is an advanced, DICOM-to-Windows printing solution with a host of rich features.
- It offers multi-modality connectivity and supports four different Windows-based printers at one time.
- For fast printing, it can be configured to send prints to the next available idle printer in rotation.
- Suitable for Pre-Natal Diagnostic Technique (PNDT) prints.
- Enhance print quality by applying linear DICOM Look-Up Tables (LUTs).
- Handles both Greyscale and colour prints at a time.
- Includes VMS OTP - An Android-based print-count recharge App for fast and timely recharge.






# VMS MediPro DICOM Software

MANAGE

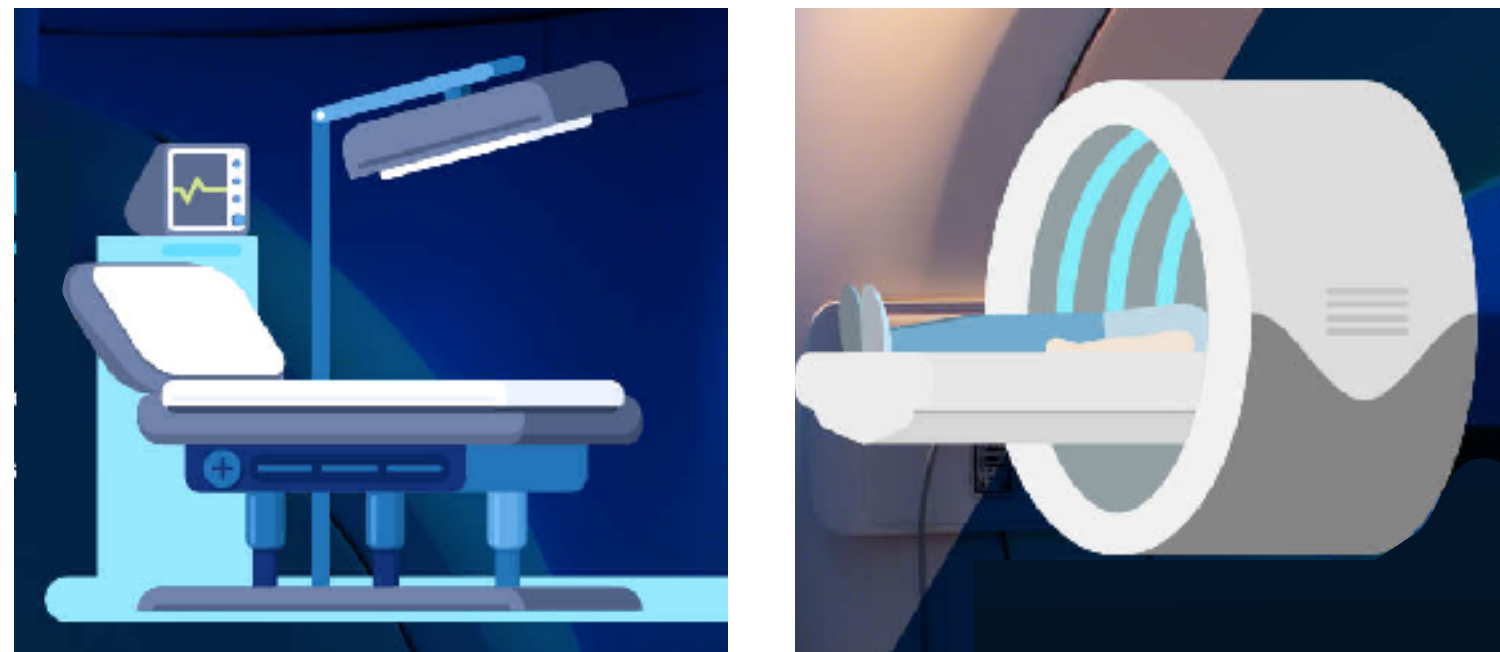

The **VMS MediPro** can receive information from unlimited different modalities

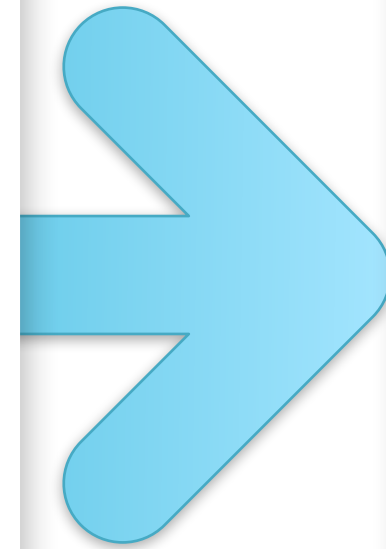

The DICOM software functions as a **receiving server** as well as a **print router** can be configured to receive information from up to **6 servers**



The received jobs can be either routed to create **PDFs** or to **print on four different inkjet / laser printers**, including in sequential printing mode




# Why Upgrade to VMS MediPro?

MANAGE

- Built on Microsoft Dot Net technology
- Supports multiple IP addresses and ports to the application, thus facilitating body-part wise LUT assignment.
- Capable of receiving both monochrome as well as colour jobs at the same time
- Superior quality prints since LUT is applied to DICOM data
- PDF is created automatically and saved in a desired folder with date and time
- Developed with licensed DICOM development tools.
- Can be upgraded to connect with new versions of modality console applications
- Get technical support for unexpected issues related to DICOM versions





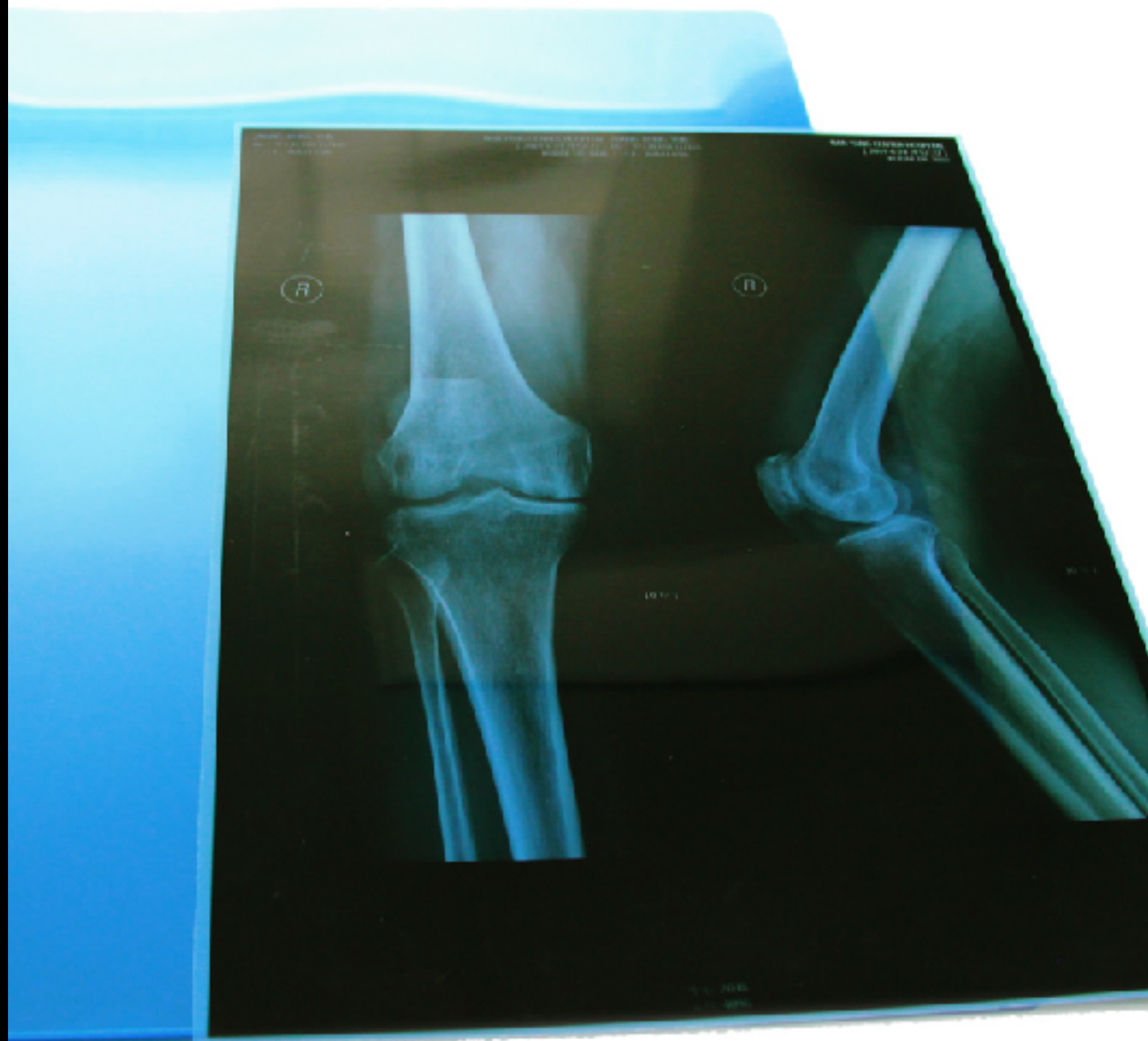

**Print**

**VMS Medical  
X-ray Dry Film**


**No environment pollution**



# VMS Medical X-Ray Dry Film

PRINT

- The **VMS Medical X-Ray Dry Film** is made from medical standard base PET material.
- Offering environmental protection coating and printable using inkjet technology, it replaces the traditional X-Ray film, saving both hospital and patients cost.
- Unique antistatic back coating to avoid sticky problem
- Water resistant and instant drying, anti-crease
- Easy storage, no need for darkroom or refrigeration, and long - term storage
- Wide and professional application in imaging department in hospital, like DR, CR, CT, MRI film output
- Available in Blue and White film base

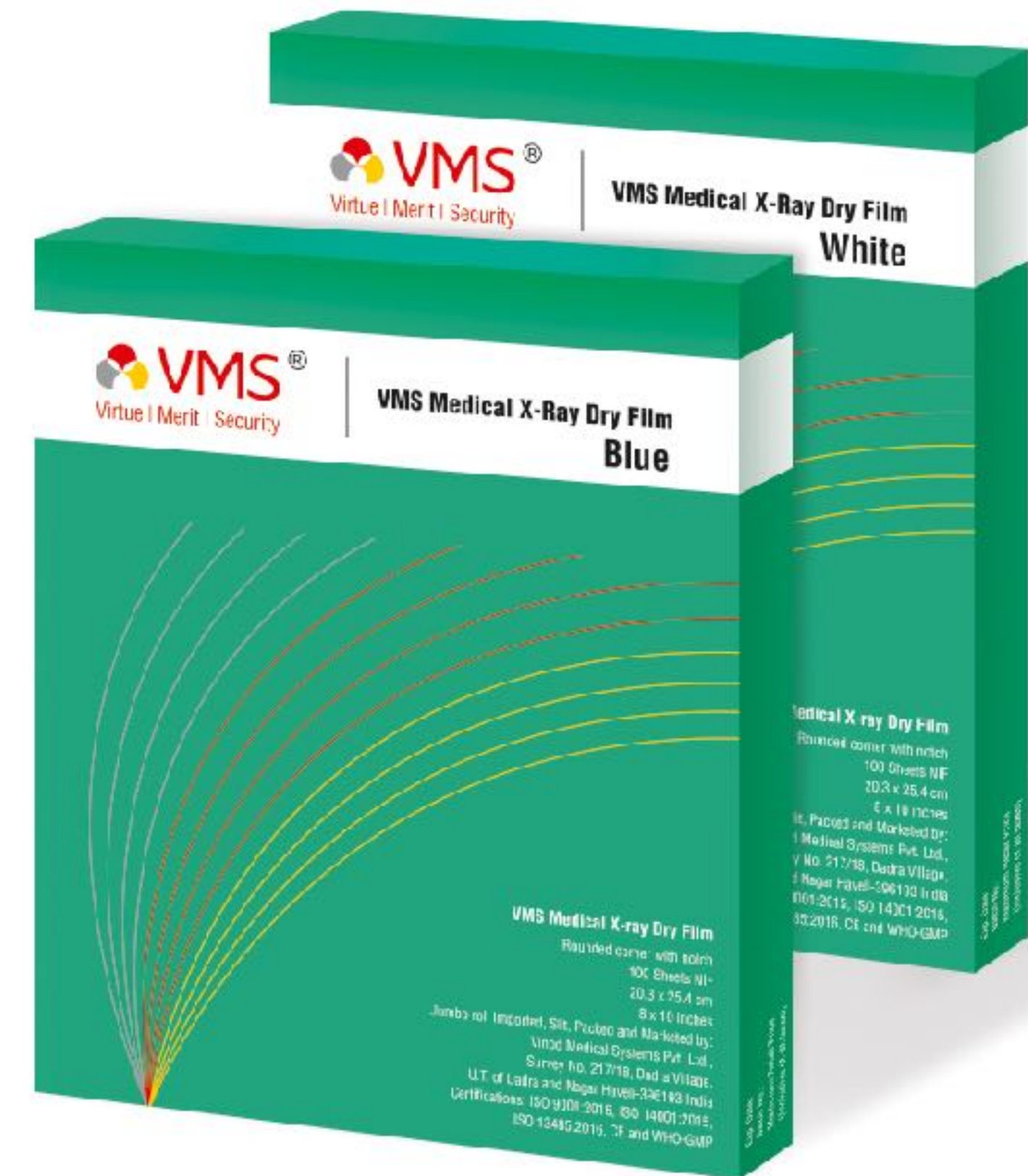





**Share**

**Instantly share  
reports with  
VMS MediPro**



# Share Reports with VMS MediPro

SHARE



- **VMS MediPro** provides automatic creation of DICOM data to hi-resolution PDFs which are saved in a desired folder with date and time
- The PDF reports can be instantly shared via WhatsApp or email
- Perfect for instant reporting requirements, especially in case of emergencies and complications
- Well-suited for sharing reports for online consultation





## CONTACT US

### CORPORATE OFFICE:

119, Omkar, The Summit Business Bay,  
1st Floor, B. L. Bajaj Road,  
Near W. E. H Metro station,  
Andheri East, Mumbai - 400 093. India.

Tel: +91 022 26820517/18

Email: [mumbai@vinodmedical.com](mailto:mumbai@vinodmedical.com)

### REGISTERED OFFICE:

Kripa Kunj, B1 - 2, Sai Nagar,  
Near Railway Crossing,  
Raipur - 492 009. India.

Tel: +91 771 4214400 / 21

Email: [vm spl@vinodmedical.com](mailto:vm spl@vinodmedical.com)



[customer care@vinodmedical.com](mailto:customer care@vinodmedical.com)

[www.vinodmedical.com](http://www.vinodmedical.com)

**Toll Free Number: 1800 572 1988**





# One Brand Many Solutions

